

10.08.2023 - Keit DLL Upgrader 1.1

This patch reduces the occurrence of software corruption during system startup and recovery, additionally reducing delay in recovery software.


Navigate to our website using the button below to download the update for installation and view the patch notes.

If you have any questions, please contact the team at support@keit.co.uk.


Please note that this patch is only for customers who have the KeitSpec IRmadillo software installed.

Steps to install the update:

1. Scroll down to "**IRmadillo™ Software Updates**" on the page linked below and click "**10.08.2023 - Keit DLL Upgrader 1.1**" to begin the download.
Once downloaded, unzip the file and navigate to the "Application" file
2. Run the file as an administrator (the patch will not install without this)

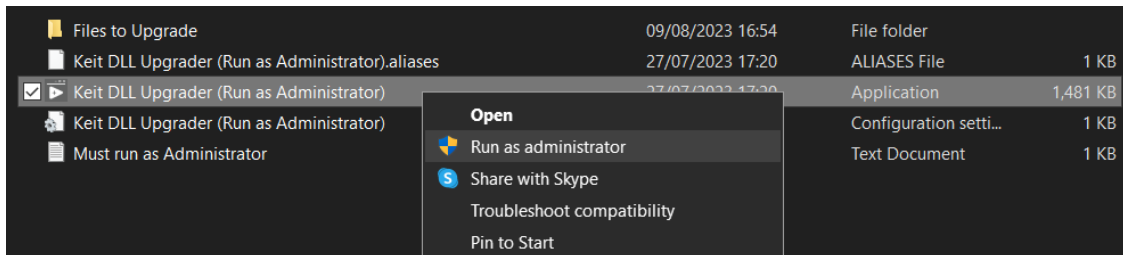

3. Enter your local admin password and the installer will automatically start updating the software. In a few moments, a message with "Upgraded Successfully". Press "Finish"

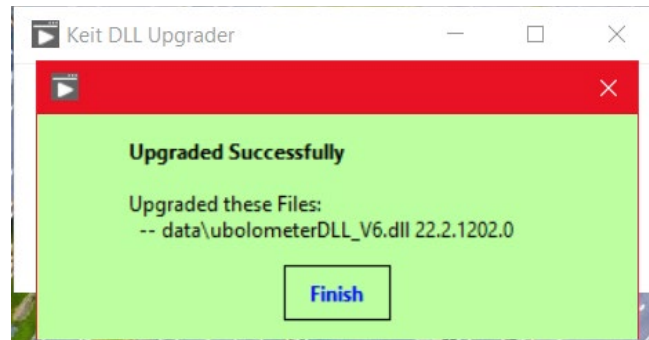

4. We recommend you reboot your computer as soon as possible after installing the update

